

# Nebulizer Instructions



1. Place compressor on a flat surface and plug into an outlet.



2. Pour medicine into nebulizer cup.



3. Screw on nebulizer cup and attach mouthpiece or mask.



4. Attach tubing to compressor and nebulizer cup.



5. Turn power switch on.



6. Inhale mist and exhale slowly.



7. Continue treatment until misting stops.



8. Turn power switch off and unplug compressor.



9. Disassemble medicine cup and rinse after each use.



10. Soak nebulizer cup and mouthpiece/mask DAILY in three parts water & one part vinegar for 30 minutes and rinse.



11. Change filter every six months.



**WRIGHT & FILIPPIS®**

(800) 482-0222

[www.FirstToServe.com](http://www.FirstToServe.com)

# Instrucciones Para El Nebulizador



1. Coloque el compresor en una superficie plana y enchufe.



2. Vierta la medicina en la taza del nebulizador.



3. Atornille la taza al nebulizador y una a la boquilla o mascara.



4. Una la tubería al compresor y a la taza del nebulizador.



5. Presione el boton de encender.



6. Inhale la niebla y exhale lentamente.



7. Siga con el tratamiento asta que la niebla pare.



8. Apague y desenchufe el compresor.




9. Desmonte la taza de medicina y enjuague despues de cada uso.



10. Una vez al dia remoje la taza y la boquilla o mascara en una mezcla de tres partes de agua y una de vinagre por 30 minutos y enjuague.



11. Cambie el filtro cada 6 meses.

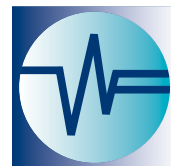

**WRIGHT  
&  
FILIPPIS®**

(800) 482-0222

[www.FirstToServe.com](http://www.FirstToServe.com)